

TERUMO PENPOL®

TRANSFUSION SYSTEMS

MULTIFUNCTIONAL HAND STRIPPER

HS100

ISO 13485 




Management System
EN ISO
13485:2012
www.tuv.com
ID 000047281

Optimize blood care with advanced technology




FEATURES

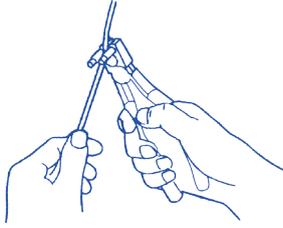


Stripping


From needle end, strip toward bag
To strip, pull tube or push rollers.

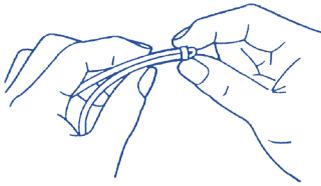

Sealing/Crimping

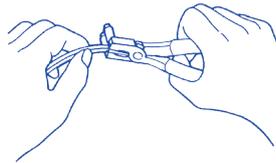

Place sealing clip onto folded tubing
between numbers, then crimp firmly
until jaws close. Check by pinching at
tubing distal end.

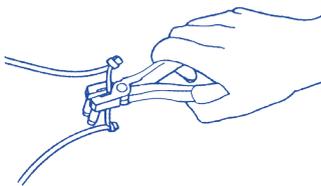

Cutting

Cut tubing between seals.

SPECIFICATIONS | MULTIFUNCTIONAL HAND STRIPPER HS100

Body metal material	Stainless Steel, 17-4PH per ASTM A747 - bead blasted & passivated
Handle grip	Plastisol, Colour-Black
Roller	Derlin AF, Colour-Black
Weight(Approx.)	Net:215g, Gross:255g
Standards & approvals	CE mark, EUroHS
Manufacturing standards	ISO 9001: 2008 , EN ISO 13485 : 2012

TERUMO PENPOL MULTIFUNCTIONAL HAND STRIPPER

- Does the work of three- strips, cuts and seals.
- Simple design for convenient usage.
- Easy handling.
- Swift and complete stripping of blood bag tubes.
- Easy-to-grip, durable handle for secure control.



REGISTERED OFFICE

TERUMO PENPOL Private Limited
I-2, JAWAHAR NAGAR
KOWDIAR P.O, TRIVANDRUM – 695003
KERALA, INDIA.
Ph.No : 0471-3015500 / 501.
Fax: 0471-2721519
CIN: U33112KL1985PTC004531

www.terumopenpol.com
www.terumobct.com
pmt@terumopenpol.com

Toll Free No : 1800 3000 7070

** Due to continuous development specifications are likely to change without prior notice
** For distribution in India only