

HS200 Multifunctional Hand Stripper ME-ACS121 Aluminum Sealing Clips



Intended Use:


The HS200 is a multifunctional tool designed to strip, crimp and cut blood bag tubes. It mechanically produces a hermetic seal with the aluminum sealing clips.

Easy Handling:

- Easy-to-grip, durable handle for secure control
- Handle design is convenient for all hand sizes


Aluminum Sealing Clips:

- Ring design for easy placement onto folded tubing
- Made of aluminum for easy, secure sealing

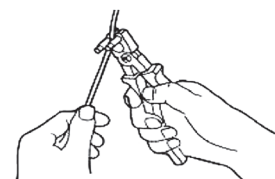


Stripping

From the needle end, strip toward the bag. To strip, pull the tubing or push the rollers.

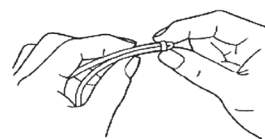


Sealing

Place the sealing clip onto the folded tubing between the numbers.

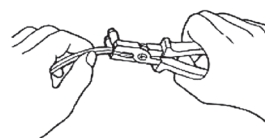

Crimping

Grasp firmly until the jaws close. Check by pinching at the distal end of the tubing.

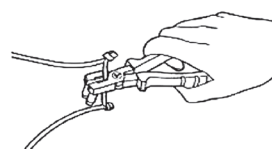

Cutting

Cut the tubing between seals.

GENERAL INFORMATION

Product codes	HS200 (Multifunctional Hand Stripper) ME-ACS121 (Aluminum Sealing Clips)
Hand stripper	1 unit per box
Sealing clips	1,000 pieces per box
Product ordering code	HS2000S000

To learn more, contact your
Terumo Blood and Cell Technologies representative.

Manufactured by Terumo Penpol Private Limited, India (hand stripper)
and Terumo Corporation, Japan (clips).
Available in select markets.

TERUMO PENPOL PRIVATE LIMITED

Bhadra Towers, TC 29/2923, Cotton Hill Road,
Vazhuthacaud, Thycaud P.O.,
Thiruvananthapuram, Kerala - 695014, India.
Phone: +91 471-3015500 / 501
E-Mail: Penpol.info@terumobct.com
TERUMOPENPOL.COM
CIN: U33112kL1985PTC004531



Terumo Blood and Cell Technologies is a medical technology company. Our products, software and services enable customers to collect and prepare blood and cells to help treat challenging diseases and conditions. Our employees around the world believe in the potential of blood and cells to do even more for patients than they do today. [TerumoBCT.com](https://www.terumobct.com)

Terumo BCT, Inc.
Lakewood, CO, USA
+1.303.231.4357

Terumo BCT Europe N.V.
Zaventem, Belgium
+32.2.715.0590

Terumo BCT Asia Pte. Ltd.
Singapore
+65.6715.3778

Terumo BCT Latin America S.A.
Buenos Aires, Argentina
+54.11.5530.5200

Terumo BCT Japan, Inc.
Tokyo, Japan
+81.3.6743.7890